

COLD


FILTER

PLUGGING

POINT

TESTER

AFP-202

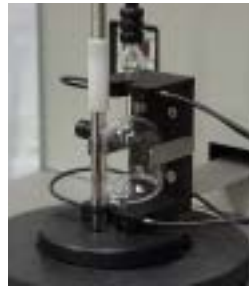

CE

with Peltier cooler

No Methanol

The model AFP-202 is a 2 test heads model of AFP series, and automatically executes the Cold Filter Plugging Point (CFPP) test method prescribed in IP309, ASTM D6371 and JIS K2288. After developing its predecessor models, this third generation model has been developed by applying Tanaka's latest technologies and the various know-hows. "Ease of operation" has been further improved on this model.

METHANOL FREE OPERATION

Use of Peltier cells for sample cooling/heating made this tester methanol free. No big chiller is needed; a small chiller with antifreeze suffices the cooling requirements.

DUAL COOLING MODES

Mode of cooling can be selected from linear cooling and conventional step-wise cooling.

SIMPLE AND EASY OPERATION

All the operation is made through a simple membrane keypad, which is standard on more than half of Tanaka's automatic petroleum testers. The operation is simple and intuitive.

RELIABLE MENISCUS DETECTION

Upper and lower meniscus detectors are mounted on a holder that slides back and forth to engage with the pipette. This new design allows not only easy operation but also reliable detection under the conditions when the pipette is frosted.

SPECIFICATIONS:

CONFORMING STANDARDS:

IP309, ASTM D6371 and JIS K2288

Type: 2 test heads.

SPECIMEN VOLUME: 45ml

MEASURING RANGE: (typical)

RT to -60 °C (The lowest temperature varies with the chiller liquid temperature and ambient temperature.)

SEQUENCE CONTROL:

By a built-in micro processor

DISPLAY:

By a fluorescence display(VFD) x 2 sets

Temperatures displayed with 0.1 °C increments

TEMPERATURE SENSOR:

PT-100

CFPP DETECTION:

By means of a photo-electric detector consisting of an LED, photo transistor, and light guide (fiber optics).

Upper and lower light guides are mounted on sliding type holder.

VACUUM REGULATOR:

Consists of a micro air pump, pressure switch, solenoid valve and a plastic bottle.

COOLING AND HEATING (needs optional chiller):

By Peltier cells. The hot side of the Peltier cells needs to be cooled by external chiller.

TEMPERATURE CONTROL:

By step(-34 °C, -51 °C and -67 °C) or linear. The linear cooling rate can be programmed.

Printer:

Built-in, thermal type.

SAFETY FEATURES:

a. Over cut: as hot side of TED reaches 60 °C , warning buzzer beeps and stops.

b. Self diagnostic: reported by an intermittent buzzer and a text message when any temperature sensor breaks down or the battery has drained out.

c. System trouble: reported by a buzzer when the built-in microprocessor runs off.

DATA OUTPUT:

RS-232C 1 channel.

POWER REQUIREMENTS:

100 or 120VAC, 10A

220 or 240VAC, 5A

Please specify AC voltage when ordering.

DIMENSIONS AND WEIGHT:

550mmWx600mmDx670mmH, 50kg

ORDERING INFORMATION:

STANDARD ACCESSORIES:

1. Temperature Sensor	2 pcs
2. Pipette	2 pcs
3. Vacuum Tubing (with Connector and Adapter)	2 pcs
4. Filter Unit (with Holder and Filter)	2 pcs
5. Filter	8 pcs
6. Lid	2 pcs
7. Test Jar	2 pcs
8. Spacers	2 pcs
9. Insulator (Large and Small)	2 sets
10. Printing Roll Paper (Thermal, W80 x D50)	1 pc
11. AC Power Cord	1 pc
12. Hose and Cramps for Chiller	1 set
13. Instruction Manual	1 copy

OPTIONAL ACCESSORIES:

Chiller for -67°C of Bath Control (ambient temp @20°C):
 TANAKA TCU-40B
 or Neslab RTE-740 (Made in USA)
 or Julabo FP40-MV (Made in Germany)

SUGGESTED SPARES FOR 2 YEARS:

1. Temperature Sensor	1 pc
2. Pipette	2 pcs
3. Vacuum Adapter (with O-Ring P-6)	1 pc
4. Vacuum Tubing (with Connector)	1 pc
5. Filter Unit (without Holder and Filter)	1 pc
6. Filter Holder (with Filter)	1 pc
7. Filter	60 pcs
8. O-Ring P-5 (for Pipette)	20 pcs
9. O-Ring P-6 (for Adapter)	10 pcs
10. O-Ring S-14 (for Filter Holder)	10 pcs
11. Test Jar	6 pcs
12. Spacers	1 pc
13. Printing Roll Paper (Thermal,W80 x D50)	20 pcs

Specifications subject to change without prior notice.

TANAKA SCIENTIFIC LIMITED

7-10-3, Ayase, Adachi-ku, Tokyo 120-0005 Japan

Tel: +81-3-3620-1711 Fax: +81-3-3620-1713

URL: <http://www.tanaka-sci.com>

e-mail: sales@tanaka-sci.com Printed in Japan 0403(E)